# RelaxBack B6 Charge

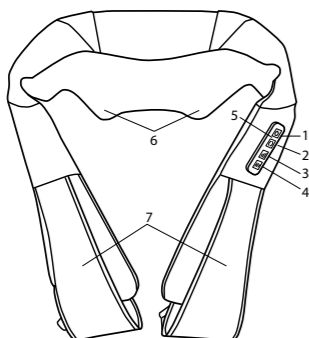
## EN – INSTRUCTIONS

### PACKAGE CONTENT

- RelaxBack B6 Charge massager
- Travel bag
- USB-C charging cable
- Manual
- Safety instructions

### PRODUCT DESCRIPTION

1. On/Off button
2. Button to change direction of rotation
3. Rotation speed button (3 speeds)
4. Heating button
5. LED charging indicator
6. Two heated rotary heads
7. Hand grips


### USE

Press the button to switch the device on. Use the button to change the direction of rotation, the button to set the 3 rotation speeds and the button to switch on heating. During use, the direction of rotation changes every minute and the device switches off automatically after 15 minutes. The device has hand grips that allow you to increase the pressure.

### Note:

- Use the device a maximum of 4 times a day. Do not use the device for more than 30 minutes at a time. Before the next use take a break of at least 30 minutes.
- In case of any health complications, consult your physician about the product's use.
- Do not use the product on areas where you feel pain or have a skin injury.
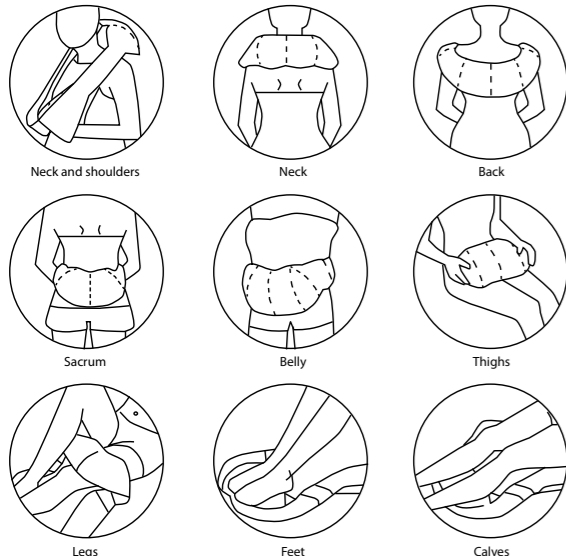
### AREAS FOR USE

### CHARGING

1. Attach the USB cable connector to the product's power cord, the other end to the adapter, and plug the adapter into a power outlet. Charge your device using the USB-C cable included.
2. The indicator lights up red during charging. Once the device is charged, the charging indicator turns green. The charging time is 4.5 hours.

*Note: For safety reasons, the device cannot be switched on and used while charging.*

### MAINTENANCE

Disconnect the device from the mains before cleaning. Do not remove the cover from the device and use a soft, damp cloth to clean it. Allow the cover to dry before use. The cover must not be washed in a washing machine.

### SPECIFICATIONS

Heating	Yes (42 °C)
Rotation in two directions	Yes
Rotation speed	3 speeds
Timer	15 minutes and automatically changes direction of rotation every minute
Number of massage heads	2 (each head has four deep-massage nodes)
Batteries	2600 mAh, Li-ion
Battery life	Up to 2 hours
Charging	USB-C cable
Charging time	4.5 hours
Input voltage	5 V – 1 A, 24 W
Dimensions	1200 x 170 x 140 mm
Weight	1477 g

### MANUFACTURER'S LIABILITY

TrueLife products have a warranty for manufacturing defects. In the event of such defects under the warranty, elem6 s.r.o. is obliged to repair or replace a defective part or product with a comparable part or product at its discretion. Rebuilding, altering or otherwise interfering with the product will void this warranty. This is the case even if the product

is lost or damaged due to a defect, negligence or other error in the product. The warranty does not cover damage resulting from use, an accident or normal wear and tear. elem6 s.r.o. assumes no liability for any accident, injury, death, loss or other claim related to or arising out of the product's use. In no event shall elem6 s.r.o. be liable for incidental or consequential damages related to or resulting from the use of the product or any of its components. Returns or exchanges of replacement parts and/or products may be subject to shipping, handling, replacement and/or restocking charges.

### SAFETY WHEN WORKING WITH THE DEVICE:

1. Read the enclosed safety information carefully.
2. For proper handling of the product, please read the current version of the manual, which is available for download on the website: [www.truelife.eu](http://www.truelife.eu).

Printing errors and changes to information reserved.

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